



Specification

Output Power	0,1 – 3,0J
Spacing	Dot mode : 1,0 – 2,0 mm (0,1 mm step) Linear mode : N/A (CW)
Length	Dot mode : 5,0 – 25 mm (5 mm step) Linear mode : 25mm
Cartridge Type	Basic (Dot/Linear) : 2,0mm / 3,0mm Essential (Dot/Linear) : 1,5mm / 2,0mm / 3,0mm Core (Dot/Linear) : 4,5mm / 6,0mm Contour (Dot/Linear) : 9,0mm / 11,0mm / 13,0mm
Dimension	1207(H) x 615,5(L) x 500(W) mm
Weight	55,1 kg



Jeisys



Adjustable Depth and Mode in a **SINGLE** cartridge

LINEAR mode

DOT mode

Design HIFU Lifting
LINEARZ

LINEARZ Only 4 Cartridges

Depth & Mode (Linear & Dot) are selectable with **just 1 cartridge**



Linear & Dot

2.0 mm
3.0 mm



Linear & Dot

1.5 mm
2.0 mm
3.0 mm



Linear & Dot

4.5 mm
6.0 mm



Linear & Dot

9.0 mm
11.0 mm
13.0 mm

LinearZ Features


Handpiece LED (White)

Mode(Linear&Dot) & Depth and current shot time can be checked in real time




Handpiece LED (BarType)

Indicate the position of the transducer and the temperature of the cartridge

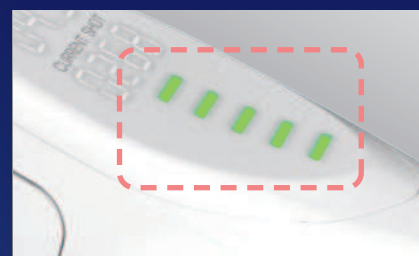

15-inch touch screen

Large display from various angles easy to secure view

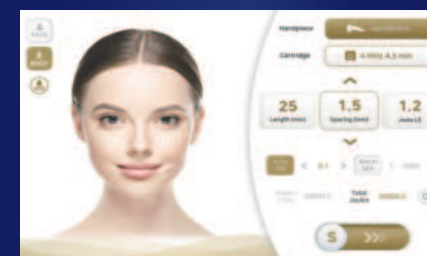

Cartridge DWR system

Equipped with a degassed water replacement system that can quickly stabilize the temperature inside the cartridge



Cartridge inner storage system

4 removable Silicone pad for cartridges
Can be stored hygienically

